


SOL-CARE™ Safety Needle

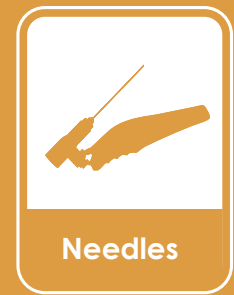

The SOL-CARE™ Safety Needle is designed to prevent needlestick injury. Its robust safety mechanism is intuitive and compatible with all standard luer lock and slip tip syringes. Available in a broad range of lengths and gauges.

- Easy, one-handed technique of activating the protective arm.
- Audible click indicates safety mechanism activation.
- Patented locking mechanism ensures needle is secure and completely locked in protective arm.
- Hub, protective arm and packaging are color coded for easy and fast identification 18G -30G needles (ISO 6009).

Two options for activation of the protective arm:



activation against hard surface




activation with thumb





SOL-CARE™ Safety Needle

- In compliance with criteria given in the **EU Directive 2010/32/EU**.
- Protective arm **doesn't obstruct visualization of injection site** (1).
- **Oriented bevel-up needle** allows for low-angle injections (2).
- Hub, protective arm and packaging are **color coded** for easy and fast identification **18G -30G needles** (ISO 6009).
- Available in 8 needle length **from 10 to 50 mm**.
- **Compatible with all standard** luer lock and slip tip **syringes**.

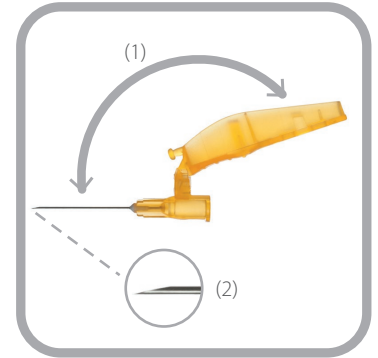

Reference number	Product description	Gauge size	Outer diameter (mm)	Length		Indications*				
				(mm)	(inch)	SC	ID	IM	IV	IA
SN3012	SOL-CARE™ Safety Needle 30G*1/2"	30G	0.3 mm	12.5 mm	1/2"	SC				
SN2715	SOL-CARE™ Safety Needle 27G*1 1/2"	27G	0.4 mm	40 mm	1 1/2"	SC	ID			
SN27125	SOL-CARE™ Safety Needle 27G*1 1/4"	27G	0.4 mm	32 mm	1 1/4"	SC	ID			
SN2734	SOL-CARE™ Safety Needle 27G*3/4"	27G	0.4 mm	18.5 mm	3/4"	SC	ID			
SN2758	SOL-CARE™ Safety Needle 27G*5/8"	27G	0.4 mm	16 mm	5/8"	SC	ID			
SN2712	SOL-CARE™ Safety Needle 27G*1/2"	27G	0.4 mm	12.5 mm	1/2"	SC	ID			
SN2612	SOL-CARE™ Safety Needle 26G*1/2"	26G	0.45 mm	12.5 mm	1/2"	SC	ID			
SN2515	SOL-CARE™ Safety Needle 25G*1 1/2"	25G	0.5 mm	40 mm	1 1/2"	SC	ID	IM		IA
SN2510	SOL-CARE™ Safety Needle 25G*1"	25G	0.5 mm	25 mm	1"	SC	ID	IM		IA
SN2558	SOL-CARE™ Safety Needle 25G*5/8"	25G	0.5 mm	16 mm	5/8"	SC	ID	IM		IA
SN2410	SOL-CARE™ Safety Needle 24G*1"	24G	0.55 mm	25 mm	1"	SC		IM		IA
SN2315	SOL-CARE™ Safety Needle 23G*1 1/2"	23G	0.6 mm	40 mm	1 1/2"			IM	IV	IA
SN23125	SOL-CARE™ Safety Needle 23G*1 1/4"	23G	0.6 mm	32 mm	1 1/4"			IM	IV	IA
SN2310	SOL-CARE™ Safety Needle 23G*1"	23G	0.6 mm	25 mm	1"			IM	IV	IA
SN2358	SOL-CARE™ Safety Needle 23G*5/8"	23G	0.6 mm	16 mm	5/8"			IM	IV	IA
SN2215	SOL-CARE™ Safety Needle 22G*1 1/2"	22G	0.7 mm	40 mm	1 1/2"			IM	IV	IA
SN22125	SOL-CARE™ Safety Needle 22G*1 1/4"	22G	0.7 mm	32 mm	1 1/4"			IM	IV	IA
SN2210	SOL-CARE™ Safety Needle 22G*1"	22G	0.7 mm	25 mm	1"			IM	IV	IA
SN2105	SOL-CARE™ Safety Needle 21G*2"	21G	0.8 mm	50 mm	2"			IM	IV	IA
SN2115	SOL-CARE™ Safety Needle 21G*1 1/2"	21G	0.8 mm	40 mm	1 1/2"			IM	IV	IA
SN21125	SOL-CARE™ Safety Needle 21G*1 1/4"	21G	0.8 mm	32 mm	1 1/4"			IM	IV	IA
SN2110	SOL-CARE™ Safety Needle 21G*1"	21G	0.8 mm	25 mm	1"			IM	IV	IA
SN2015	SOL-CARE™ Safety Needle 20G*1 1/2"	20G	0.9 mm	40 mm	1 1/2"			IM	IV	
SN2010	SOL-CARE™ Safety Needle 20G*1"	20G	0.9 mm	25 mm	1"			IM	IV	
SN1920	SOL-CARE™ Safety Needle 19G*2"	19G	1.1 mm	50 mm	2"			IM	IV	
SN1915	SOL-CARE™ Safety Needle 19G*1 1/2"	19G	1.1 mm	40 mm	1 1/2"			IM	IV	
SN1910	SOL-CARE™ Safety Needle 19G*1"	19G	1.1 mm	25 mm	1"			IM	IV	
SN1820	SOL-CARE™ Safety Needle 18G*2"	18G	1.2 mm	50 mm	2"				IV	
SN1815	SOL-CARE™ Safety Needle 18G*1 1/2"	18G	1.2 mm	40 mm	1 1/2"				IV	
SN1810	SOL-CARE™ Safety Needle 18G*1"	18G	1.2 mm	25 mm	1"				IV	

SC - subcutaneous, ID - intradermal, IM - intramuscular, IV - intravenous, IA - intra-arterial.

* **Recommended indications for specific measures. Users are explicitly responsible for correct application.**

Sales Unit: 100 pcs/Box; 800 pcs/Case

For more information, please contact:

Sol-Millennium Europe Sp. z o.o.
 Twarda 18, PL-00105 Warsaw, Poland
 Tel: +48 22 202 67 23, e-mail: info.eu@sol-m.com

sol-m.com

2018-03-SN-0224

