





FX Wire® Advanced Measurement Guidewire

Designed to provide easy access to challenging biliary anatomy, these guidewires combine a clear measurement system with a flexible, 5cm hydrophilic coated radiopaque distal tip for a simple and streamlined clinical experience.

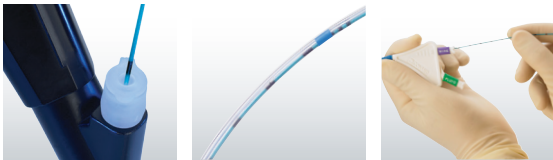

APOLLO 3 AC™
TRIPLE LUMEN PAPILOTOMES

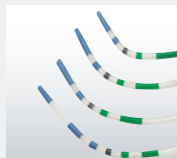

PROFORMA™
HF 4.5 CANNULA




DURAGLIDE 3™
TRIPLE LUMEN BALLOON



To learn more about these and other innovative products, call **800-448-6506** or visit **ConMed.com**.

ADVANCED
MEASUREMENT
GUIDEWIRE

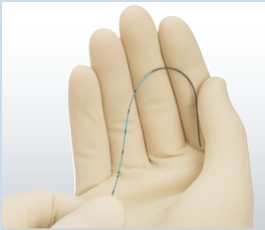

KINK RESISTANT

Benefits Beyond Measurement

When it comes to device delivery in tough-to-reach biliary anatomy, our FX Wire® fully-integrated measurement guide-wires provide the size and simplicity necessary for easy instrument access and guidance.

Hydrophilic, Flexible, Atraumatic Tip

Facilitates access past obstructions and tortuous anatomy

REGILIAN™ Nitinol Core

Kink-resistant, optimal column strength provides enhanced stricture access, easy device exchange and superior torque response


Lubricious Silicone Coating

Provides smooth, seamless device exchange

Light Blue PTFE Jacket with Black Markings

Enhances visualization by reducing endoscopic glare

FX WIRE® FEATURES

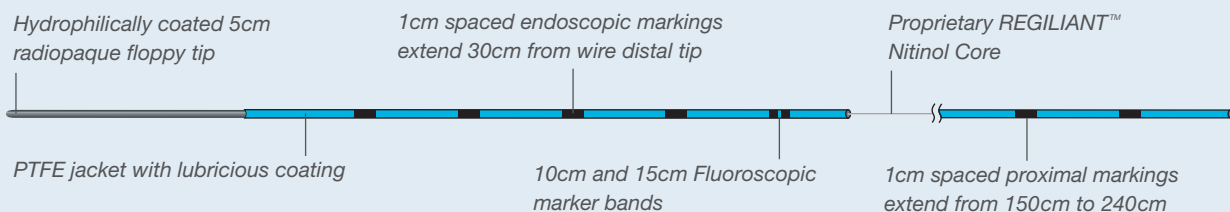

Endoscopic Measurement Made Easy

Endoscopic Measurement System


Individual 1cm spaced markings and multiple markings every 5cm for ease of measurement and to help detect inadvertent movement

Radiopaque Marker Bands

Featuring Radiopaque Marker Bands at 10cm and 15cm, which helps facilitate accurate fluoroscopic stricture measurement

To order or learn more, contact
your local ConMed Representative
or visit **ConMed.com**

ADVANCED
MEASUREMENT
GUIDEWIRE

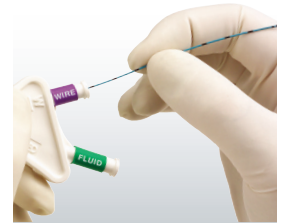

FULLY INTEGRATED
MEASUREMENT
GUIDEWIRE



525 French Road, Utica, NY 13502

ADVANCED
MEASUREMENT
GUIDEWIRE

Ordering Information

Description	Guidewire O.D.	Length	Tip Type	Torque Vise	Packaging	Catalog Number
Standard Wire	.035" / .89 mm	450 cm	Straight Tip	No	3/box	FX35450S
Standard Wire	.035" / .89 mm	450 cm	Torquable Angled Tip	Yes	3/box	FX35450A
Standard Wire	.035" / .89 mm	260 cm	Straight Tip	No	3/box	FX35260S
Standard Wire	.035" / .89 mm	260 cm	Torquable Angled Tip	Yes	3/box	FX35260A

Customer Service: 1-800-448-6506
International Sales: 1-315-797-8375
Fax: 1-800-438-3051

ConMed.com
info@mail.conmed.com